


**NOTES:**

Opening suits a single SlidaScope Cabinet with IAS Fitted. Power; Standard 240V 10A GPO/DGPO (Recommended) 'Type I' DGPO's Shown.

One each of DGPO & Data Sockets.

50mm gap between opening & nearest part of DGPO/Data Plates to allow room for Flashings supplied with Cabinet.

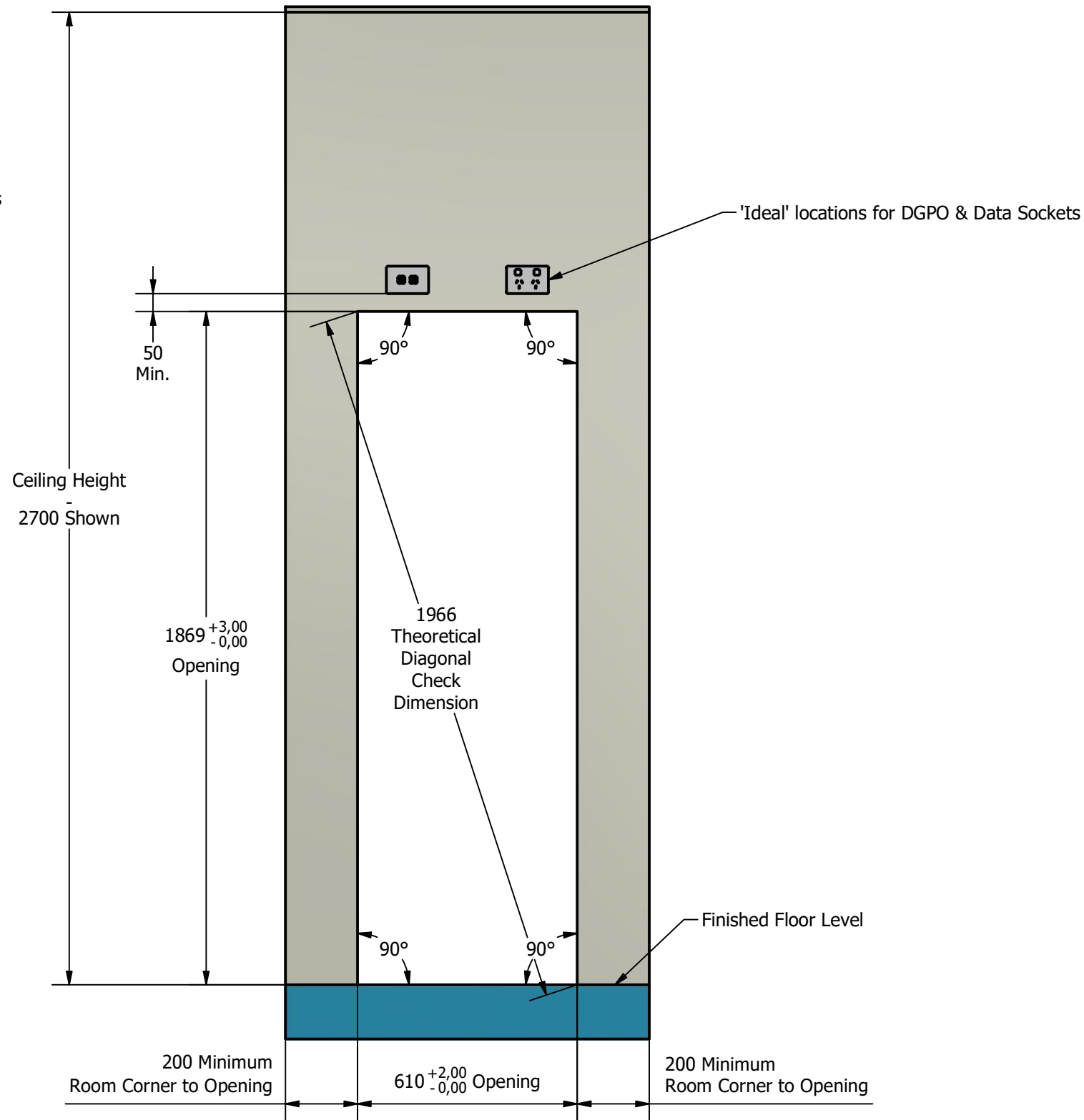

**TOLERANCE NOTE**

**LINEAR DIMENSIONS;**  
 WHOLE DIMENSION = ±0.25mm  
 DIMENSION 0.0 = ±0.20mm  
 DIMENSION 0.00 = ±0.10mm

**ANGULAR DIMENSIONS;**  
 WHOLE DIMENSION = ±0.5°  
 ANGULAR 0.0° = ±0.25°

**UNLESS NOTED OTHERWISE**


FINISHING INSTRUCTIONS: <b>ALL EDGES &amp; HOLES DEBURRED - NO RAW EDGES WILL BE ACCEPTED</b>		55 Cordwell Rd YANDINA QLD 4561 <b>Tel:</b> 07 5478 9977 <b>Fax:</b> 07 5478 8854 <b>Website:</b> www.smartlinemedical.com		AUTHOR: Frank Grabs	DESCRIPTION: Wall Opening, Power & Data Locations - Single Pass-Through SlidaScope Cabinet	TITLE: SlidaScope Cabinets
WORK TO BE COMPLETED BY:	TOLERANCE NOTE:	 ABN 34 119 209 611		CHECKED:	DXF FILES:	SALES ORDER / JOB ORDER / OSTENDO BOM CODE:
REVISION REQUEST:	MATERIAL:			SHEET: 1 OF 1		PART NUMBER: <b>WallOpening Slida1</b>
<b>A3</b>	LABOUR QUOTED/COMPLETED:			REVISION: A	DATE: 15/09/2021	
	© Copyright of Smartline Medical P/L Unauthorised Duplication or Transmission in any Form is illegal & may result in Legal Action					

**NOTES:**

Opening suits Two SlidaScope Cabinets with IAS Fitted.  
Power; Standard 240V 10A GPO/DGPO (Recommended)  
'Type I' DGPO's Shown.

Two each of DGPO & Data Sockets.

50mm gap between opening & nearest part of DGPO/Data Plates  
to allow room for Flashings supplied with Cabinet.

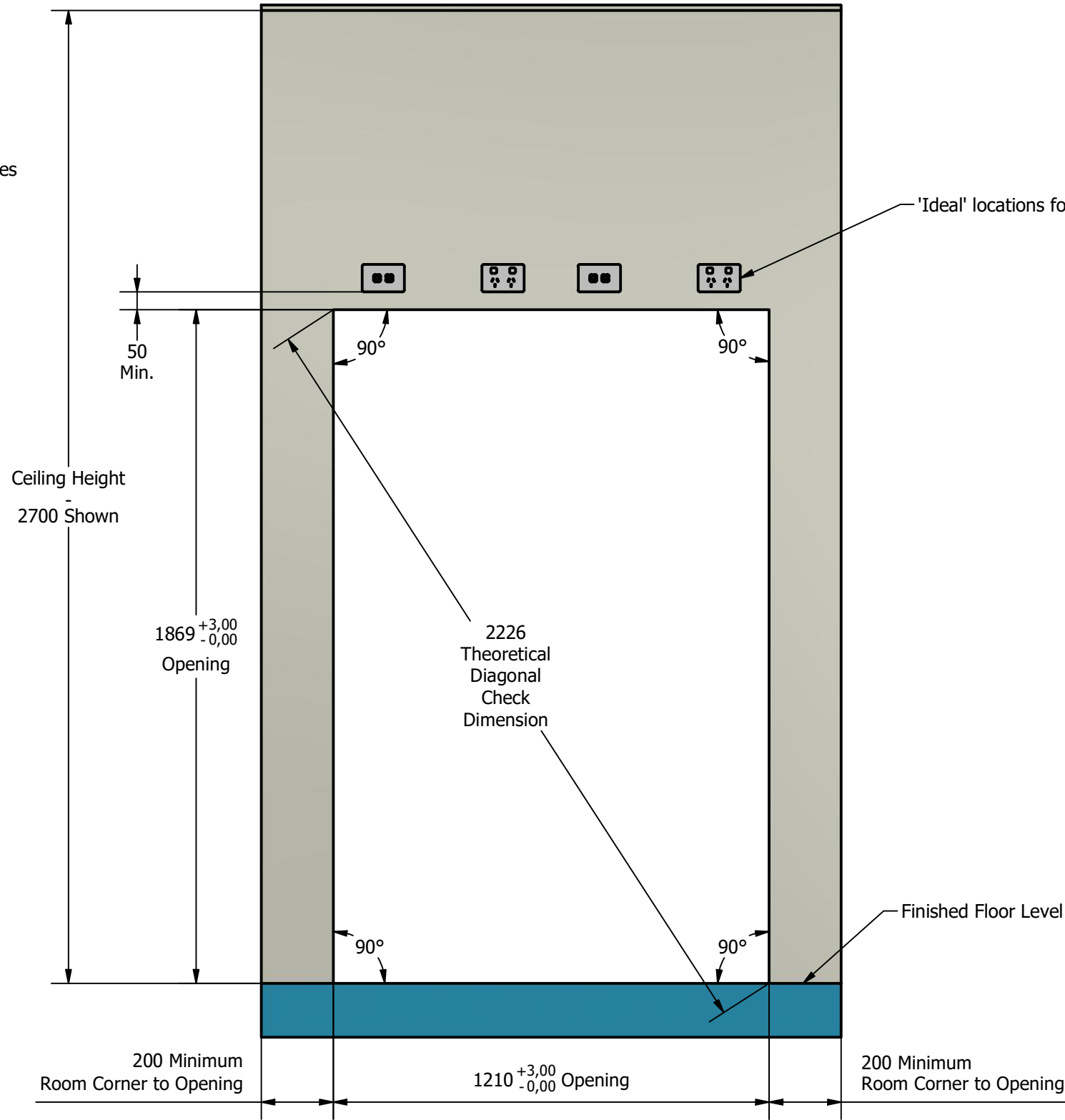

**TOLERANCE NOTE**

**LINEAR DIMENSIONS;**  
WHOLE DIMENSION = ±0.25mm  
DIMENSION 0.0 = ±0.20mm  
DIMENSION 0.00 = ±0.10mm

**ANGULAR DIMENSIONS;**  
WHOLE DIMENSION = ±0.5°  
ANGULAR 0.0° = ±0.25°

**UNLESS NOTED OTHERWISE**

FINISHING INSTRUCTIONS: <b>ALL EDGES &amp; HOLES DEBURRED - NO RAW EDGES WILL BE ACCEPTED</b>		55 Cordwell Rd YANDINA QLD 4561 <b>Tel:</b> 07 5478 9977 <b>Fax:</b> 07 5478 8854 <b>Website:</b> www.smartlinemedical.com		AUTHOR: Frank Grabs	DESCRIPTION: Wall Opening, Power & Data Locations - Twin Pass-Through SlidaScope Cabinets	TITLE: SlidaScope Cabinets
WORK TO BE COMPLETED BY:	TOLERANCE NOTE:	 ABN 34 119 209 611		CHECKED:	DXF FILES:	SALES ORDER / JOB ORDER / OSTENDO BOM CODE:
REVISION REQUEST:	MATERIAL:			SHEET: 1 OF 1		PART NUMBER: <b>WallOpening Slida2</b>
<b>A3</b>	LABOUR QUOTED/COMPLETED:			REVISION: A	DATE: 15/09/2021	


© Copyright of Smartline Medical P/L  
Unauthorised Duplication or Transmission in any Form  
is illegal & may result in Legal Action

**NOTES:**

Opening suits Three SlidaScope Cabinets with IAS Fitted. Power; Standard 240V 10A GPO/DGPO (Recommended) 'Type I' DGPO's Shown.

Three each of DGPO & Data Sockets.

50mm gap between opening & nearest part of DGPO/Data Plates to allow room for Flashings supplied with Cabinet.

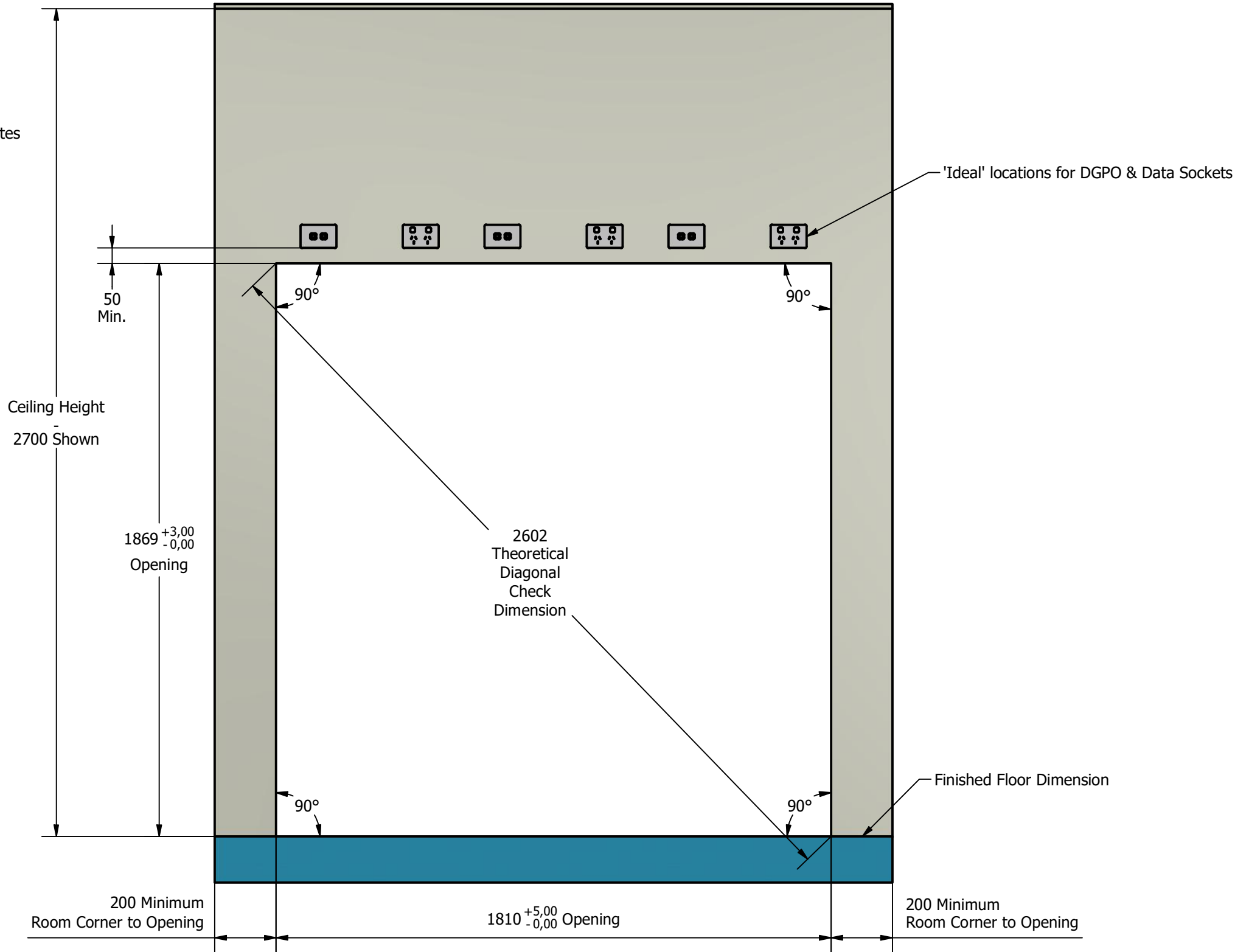

**TOLERANCE NOTE**

**LINEAR DIMENSIONS;**  
 WHOLE DIMENSION = ±0.25mm  
 DIMENSION 0.0 = ±0.20mm  
 DIMENSION 0.00 = ±0.10mm

**ANGULAR DIMENSIONS;**  
 WHOLE DIMENSION = ±0.5°  
 ANGULAR 0.0° = ±0.25°

**UNLESS NOTED OTHERWISE**

FINISHING INSTRUCTIONS: <b>ALL EDGES &amp; HOLES DEBURRED - NO RAW EDGES WILL BE ACCEPTED</b>		55 Cordwell Rd YANDINA QLD 4561 <b>Tel:</b> 07 5478 9977 <b>Fax:</b> 07 5478 8854 <b>Website:</b> www.smartlinemedical.com		AUTHOR: Frank Grabs	DESCRIPTION: Wall Opening, Power & Data Locations - Triple Pass-Through SlidaScope Cabinets	TITLE: SlidaScope Cabinets
WORK TO BE COMPLETED BY:	TOLERANCE NOTE:	 ABN 34 119 209 611		CHECKED:	DXF FILES:	SALES ORDER / JOB ORDER / OSTENDO BOM CODE:
REVISION REQUEST:	MATERIAL:			SHEET: 1 OF 1		PART NUMBER: <b>WallOpening Slida3</b>
<b>A3</b>	LABOUR QUOTED/COMPLETED:			REVISION: A	DATE: 15/09/2021	
	© Copyright of Smartline Medical P/L Unauthorised Duplication or Transmission in any Form is illegal & may result in Legal Action					