

**SMARTLINE**  
MEDICAL

# TOE Probe Cabinet

CLASSIC SERIES

---



# SERIES OPTIONS

## Core Series



**RotaScope Core**  
ROTA-CO-T

## CLASSIC SERIES

---



**RotaScope Classic**  
ROTA-T-S

# OVERVIEW

Feature	Core	Classic
<b>Integrated Air System</b> Fully independent air system. IAS requires no infrastructure from the hospital and provides clean HEPA filtered air to the cabinet.	✓	✓
<b>Plug and Play</b> No complex install requirements	✓	✓
<b>Customisable Configuration</b> (single door / Pass-Thru)	✓	✓
<b>Low Footprint – High Density Storage</b> Stores 9 probes in just 60cm x 60cm footprint		
<b>EN16442 Compliant</b>		✓
<b>Probe Data Monitoring &amp; Intelligent Environmental Monitoring</b> Touch Screen Display shows remaining probe storage time and environmental monitoring		✓
<b>Intelligent Interior LED Lighting</b>		✓
<b>Remote Log In</b> 24/7 service capability and reporting		✓



# SECURE HANGING

- Top hangars provide secure probe storage
- Gates to prevent Probe from falling out
- Half moon on top hangar to prevent crimping of probe




# SECURE HANGING



- Bottom cup mount designed to hold probe video blocks of all types and sizes
- Cable management slots prevent umbilical cord touching insertion tube
- Bottom mounts can be adjusted in height for shorter probes such as paediatric

# INTEGRATED AIR SYSTEM

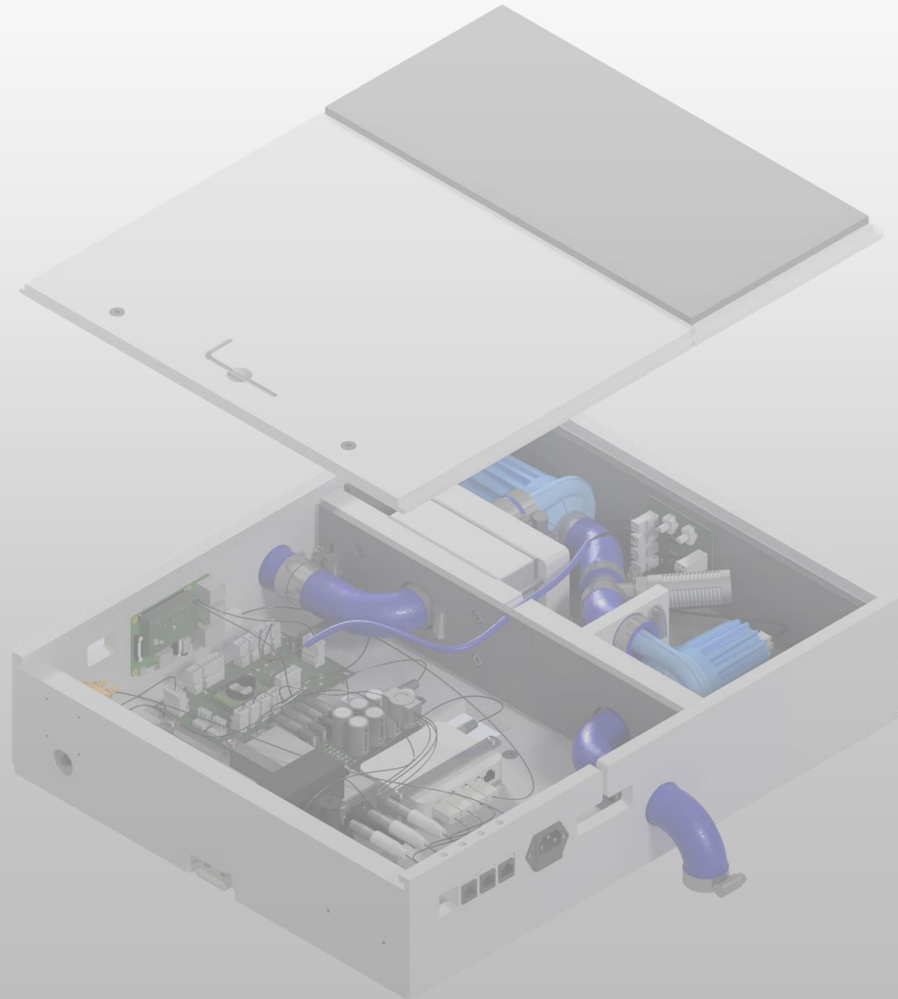

- Fully independent air system
- Requires no infrastructure from the hospital
- Quiet
- Service required only once annually

# CLASSIC SERIES

## TRACKING AND TRACEABILITY



- Barcode tracking system
- Touch screen display shows remaining probe validity storage time in and out of cabinet
- Compatible with any washing method
- Paperless validation of process

The Environmental Monitor measures and records the following cabinet conditions:

- Temperature
- Humidity
- Air (Cabinet Pressure)
- Door Status
- IAS Temperature
- Routine Cleaning Status

# CLASSIC SERIES

## TICKET PRINTOUT

### SCOPE PROCESSING & STORAGE RECORD

SCOPE TYPE: ENDOSCOPE  
SCOPE NAME: MAJ-1193  
WASHER: AER 2  
WASHED AT: 14:15:21 Nov-21 18  
REMOVED AT: 14:18:19 Nov-21 18  
STORED IN: ROTASCOPE A  
OPERATOR: door user  
SINCE WASH: 0.1 HOURS  
OOC EXPIRES: 17:16:08 Nov-21 18

- Probes can be scanned at printer at any stage for instant status report
- Tickets can be customised based on site preference
- Printed onto stickers for easy physical record keeping and daybook procedures

### SMART STATION

Plus 2 Cables  
One Station  
does everything with 10m Radius





**SMARTLINE**  
MEDICAL

# TOE Probe Cabinet

CLASSIC SERIES

---

